SIMPLIFIED USER GUIDE VITEK 2

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
Card Type:

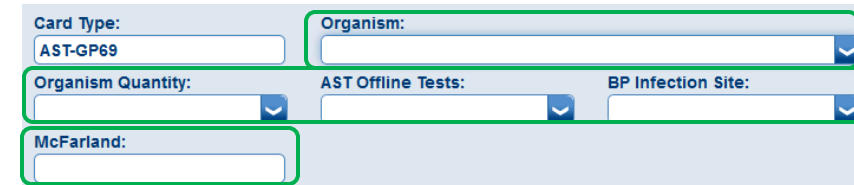Card Type: Organism:

Organism Quantity: AST Offline Tests: BP Infection Site:

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The screenshot displays the VITEK FLEXPREP software interface. At the top, there are two main tabs: "Sample Preparation" (yellow) and "Configuration" (dark blue). The "Configuration" tab is active. The interface includes a header with the BIOMÉRIEUX logo. Below the header, there are input fields for "Lab ID:" and "Isolate:". The main section is titled "VITEK® 2" and contains a "VITEK 2 PC:" label and a "Cassette ID:" label. The "Cassette ID:" field is highlighted with a red rectangle and contains the text "URI". To the right of the "Cassette ID:" field, there is a button labeled "Quick McFarland Entry". Below the "Cassette ID:" field, there is a section titled "Instruments" which contains two identical entries for "Simulator" with "Available Slots: 120". On the right side of the interface, there are three buttons: "Layout", "Validate" (with a checkmark icon), and "Cancel" (with an 'X' icon). The bottom status bar shows the user "labsuper", a help icon, and the text "VITEK FLEXPREP Release version 1" and "14:10".

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The screenshot displays the VITEK FLEXPREP software interface. At the top, there are two main tabs: "Sample Preparation" (yellow) and "Configuration" (dark blue). The "Sample Preparation" tab is active. Below the tabs, the "BIOMÉRIEUX" logo is visible on the right. The main workspace is divided into several sections. On the left, there is a "VITEK® 2" section with a "Summary" tab. Below this, a table lists the setup for the VITEK 2, with columns for Slot #, QC, Accession ID, Card, and Organism. The first row (Slot 1) is highlighted in yellow and shows "103-1" in the Accession ID column. To the right of the table, there are input fields for "Lab ID:" (containing "103") and "Isolate:" (containing "1"). Further right, there are fields for "Patient:", "Card Type:", "Organism:", "Organism Quantity:", "AST Offline Tests:", "BP Infection Site:", and "McFarland:". Below these fields, there is a visual representation of a VITEK 2 cassette with a yellow circle labeled "1" indicating the first slot. At the bottom of the main workspace, there is a "Cassette ID:" field with the value "URI". On the right side of the interface, there is a vertical toolbar with buttons for "Layout", "New Cassette" (F6), "Add Card" (F8), "Delete Card", "McFarland Ent..." (F12), "Summary" (F9), and "Send Cassette" (F10). At the bottom of the screen, there is a status bar with the text "VITEK FLEXPREP Release version 1" and the time "14:13".

Slot #	QC	Accession ID	Card	Organism
1		103-1		
2				
3				
4				
5				
6				
7				
8				
9				
10				

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Patient: Schleicher, Butch

Card Type: AST-GN04 Organism: Escherichia coli

Organism Quantity: moderate (mod2) AST Offline Tests: BP Infection Site: Urine

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Senden des Carriers


The screenshot displays the VITEK FLEXPREP software interface. At the top, there are tabs for 'Sample Preparation' and 'Configuration', and the 'BIOMÉRIEUX' logo. Below these, there are input fields for 'Lab ID:' and 'Isolate:'. The main area is divided into several sections:

- VITEK® 2 Summary:** A table with columns 'Slot #', 'QC', 'Accession ID', 'Card', and 'Organism'. The table contains two rows of data, with the third row highlighted in yellow.
- QC:** A section for Quality Control with a 'Patient:' label and various input fields for 'Card Type:', 'Organism:', 'Organism Quantity:', 'AST Offline Tests:', 'BP Infection Site:', and 'McFarland:'.
- PC: VITEK:** A section with buttons for 'New Cassette' (F6), 'Add Card' (F8), 'Delete Card', 'McFarland Ent...', 'Summary' (F9), and 'Send Cassette' (F10). The 'Send Cassette' button is highlighted with a red box.
- Layout:** A button in the top right corner.
- Validate:** A button with a checkmark icon.
- Cancel:** A button with an 'X' icon.

At the bottom, there is a status bar with the text 'VITEK FLEXPREP Release version 1' and the time '14:23'.

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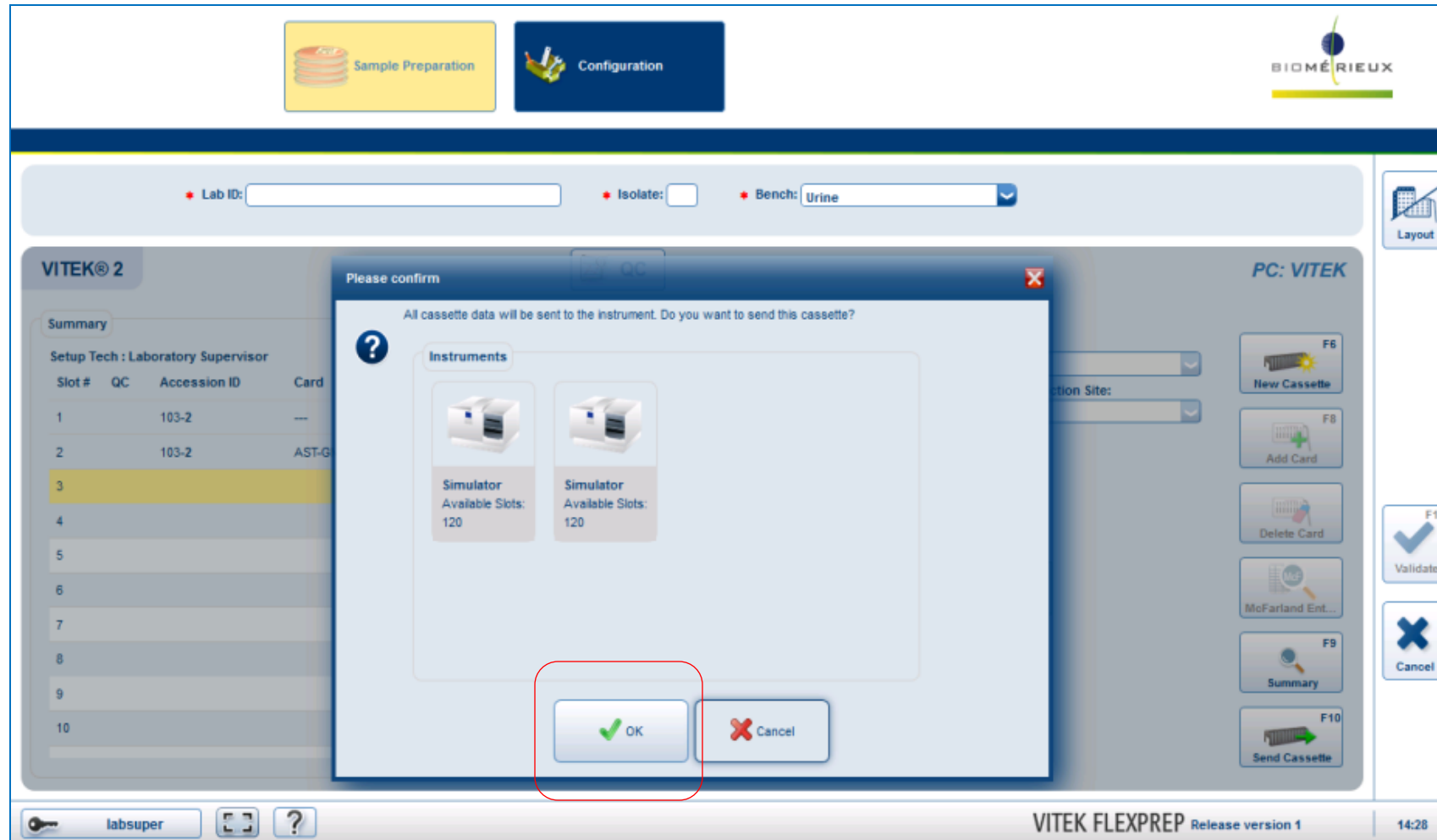The screenshot displays the VITEK FLEXPREP software interface. At the top, there are two main tabs: "Sample Preparation" (yellow) and "Configuration" (dark blue). The "Configuration" tab is active. Below the tabs, there are input fields for "Lab ID:", "Isolate:", and "Bench: Urine". The main area is divided into a "Summary" section on the left and a "Patient" section on the right. The "Summary" section contains a table with columns: Slot #, QC, Accession ID, Card, and Organism. The table has 10 rows, with row 3 highlighted. The "Patient" section contains fields for "Card Type:", "Organism:", "Organism Quantity:", "AST Offline Tests:", and "BP Infection Site:". A central dialog box titled "Information" is open, displaying the message "Cassette information successfully sent to VITEK 2 Systems." with a green checkmark and an "OK" button. On the right side of the interface, there are several function buttons: "New Cassette" (F6), "Add Card" (F8), "Delete Card" (F12), "Validate" (checkmark), "Cancel" (X), "Summary" (F9), and "Send Cassette" (F10). The bottom status bar shows the user "labsuper", the version "VITEK FLEXPREP Release version 1", and the time "14:34".

Slot #	QC	Accession ID	Card	Organism
1		103-2	---	
2		103-2	AST-GN04	
3				
4				
5				
6				
7				
8				
9				
10				

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Organism: Esch.coli

AES Findings: Consistent

Phenotypes Selected for Review: AMINOGLYCOSIDES, RESISTANT GEN TOB NET AMI (AAC(5)?), BETA-LACTAMS, EXTENDED SPECTRUM BETA-LACTAMASE

AST-V102		GN	
2 - APPA	3 - ADO	12 - ACITp	13 - dGLU
10 - H2S	11 - RNAG	19 - dMAN	20 - dVME
17 - BGLU	18 - dMAL	27 - PLE	29 - TyA
23 - PwA	26 - LP	35 - dTRE	36 - CIT
33 - SAG	34 - dTAG	42 - SUCT	43 - NAGA
40 - ILATk	41 - AGLU	48 - LDC	53 - IHISa
46 - GlyA	47 - DDC	61 - IMLi	62 - CLLM
50 - O129R	59 - GSA		64 - LATa

Patient Name: Wilson, Kelsey

Analysis Status: 5.00 hr - Final

Required ID Online tests: None required

Analysis Messages: ID Confidence: Excellent identification

Ulnumber: 0401411040524500

McFarland: []

Setup Tech: Laboratory Supervisor (Lab Su...)

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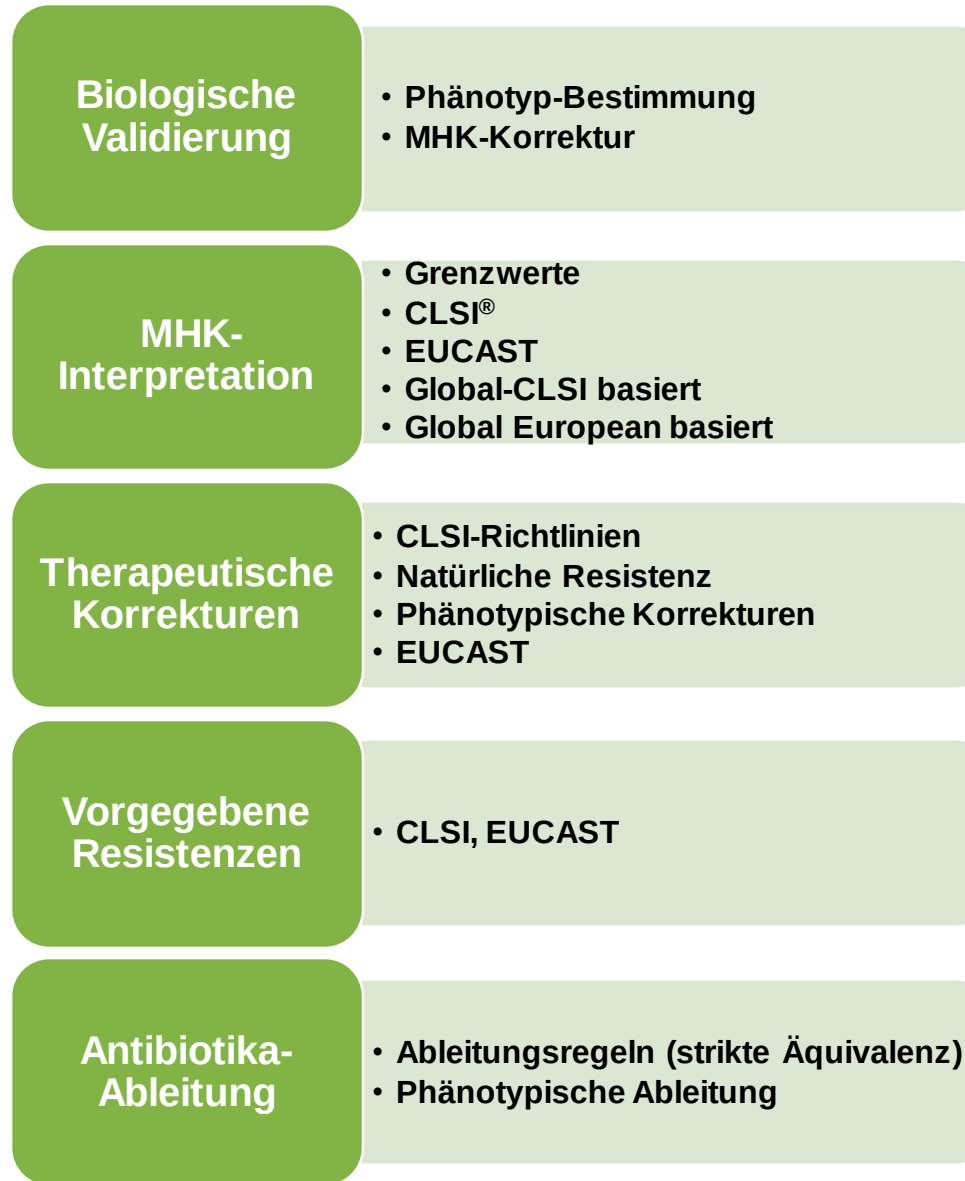
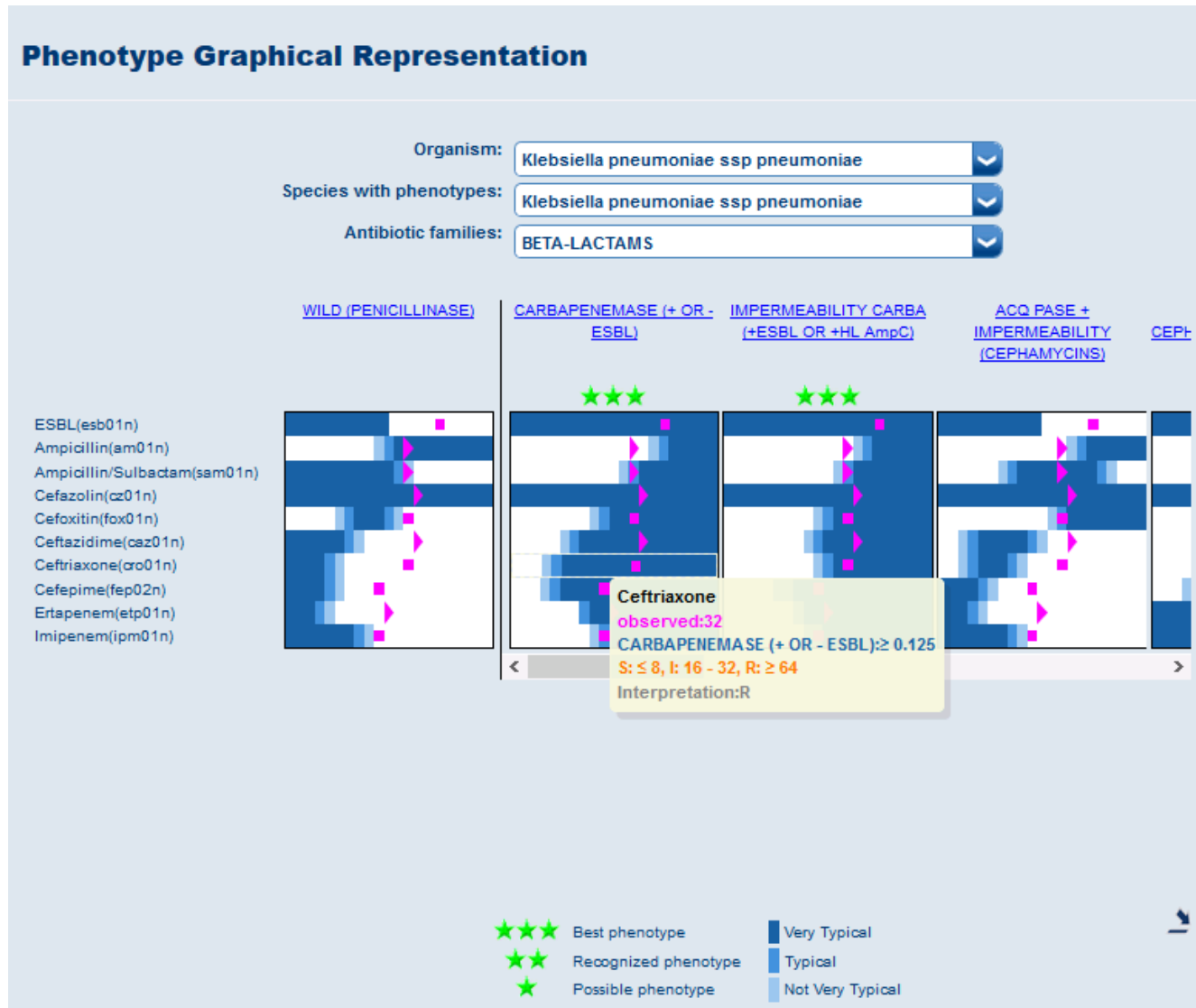
Shipments

Lot Number: 1491114403 Expiration Date: 12/14/2020
Date Received: 5/26/2019 Card Type: AST-XN06
Quantity: 1
Certificate of Conformance: ☒

Add Shipment

Date Received	Lot Expiration	Quantity	Lot Number	Card Type	Certificate of Conformance	Receiver ID	Edit
5/29/2019	5/29/2020	4	2410914203	QN	☒	Labadrin	
5/29/2019	12/22/2020	1	6851122103	AST-GH95	☒	Labadrin	
5/17/2019	12/14/2020	14	2431104103	YET	☒	Labadrin	

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VITEK® 2 **QC** **PC: VITEK**

Summary

Setup Tech : Laboratory Administrator

Slot #	QC	Accession ID	Card	Organism
1	✓		GN	
2				
3				
4				
5				
6				
7				
8				
9				
10				

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ANNEX 2: EINGABE NEUER AST-KARTEN

Maintain AST Card Definitions

	Card Type	Card Class
☐	AST-EKX3	Gram Negative Susceptibility
☐	AST-EKX4	Gram Negative Susceptibility
☐	AST-EKX5	Gram Negative Susceptibility
☐	AST-EKX6	Gram Negative Susceptibility

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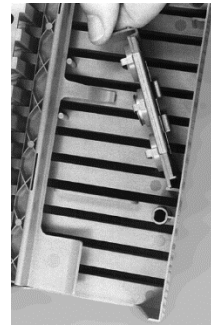
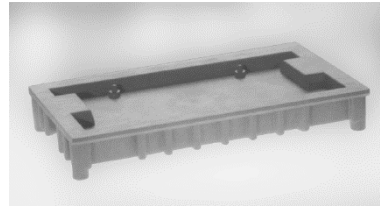
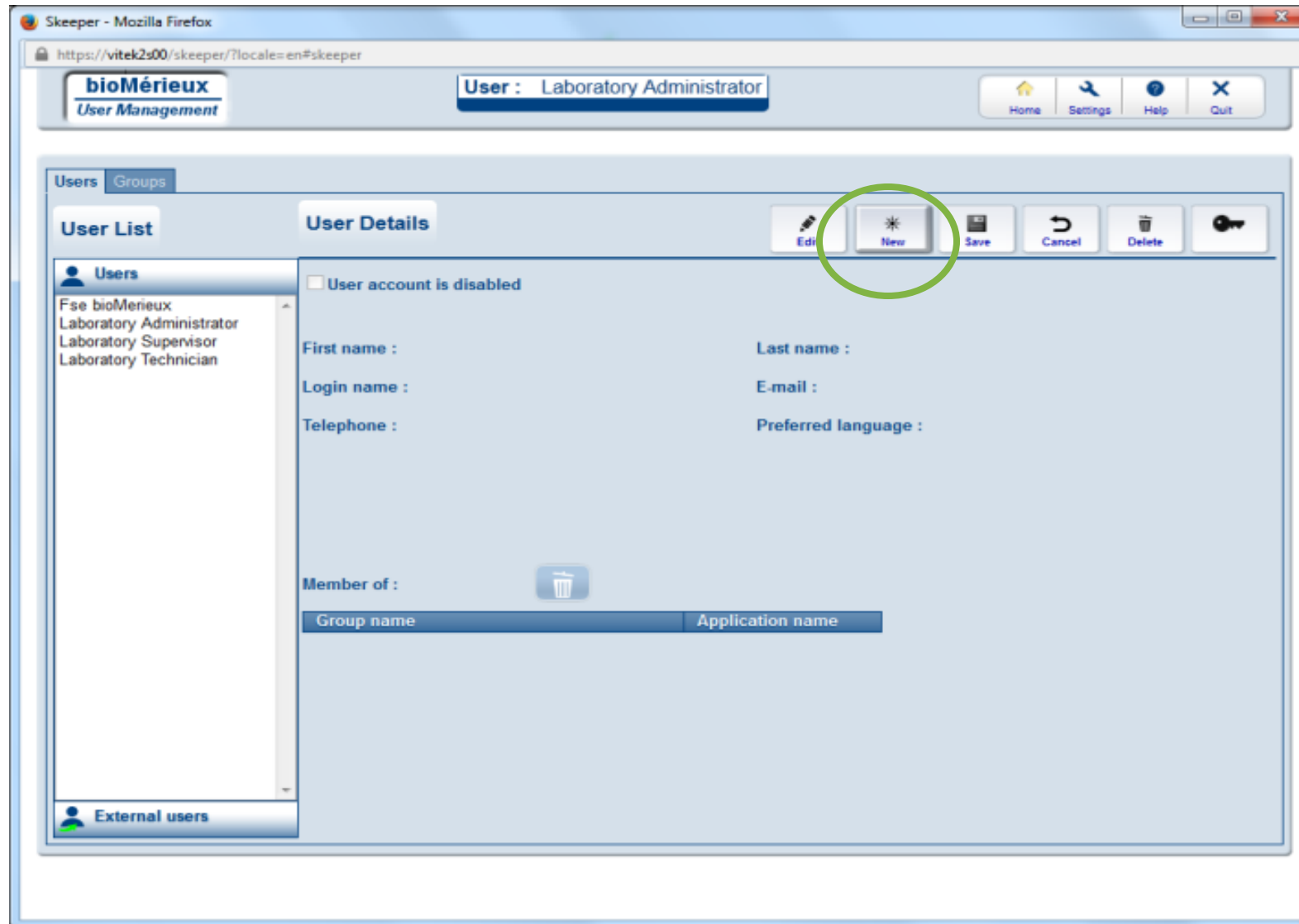
User List

Users

Fse bioMérieux
Laboratory Administrator
Laboratory Supervisor
Laboratory Technician

User Details

☐ User account is disabled

First name : [red field] Last name : [red field]

Login name : [red field] E-mail : [white field]

Telephone : [white field] Preferred language : [dropdown menu]

Password : [red field] Confirm password : [red field]

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Screenshot of the bioMérieux User Management interface in a Mozilla Firefox browser. The interface shows the 'User Details' tab for a user named 'John Doe'. The 'Preferred language' dropdown menu is open, displaying a list of languages including Deutsch, Ελληνικά, English, español, français, italiano, 日本語, 한국어, polski, português, português do Brasil, русский, Türkçe, and 中文. The 'User List' on the left shows a list of users: Fse bioMérieux, Laboratory Administrator, Laboratory Supervisor, Laboratory Technician, and John Doe. The 'Member of' section shows a trash icon. The 'External users' section is at the bottom left.

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
Speichern klicken

The screenshot shows the 'bioMérieux User Management' web application. The browser title is 'Skeeper - Mozilla Firefox' and the address bar shows 'https://vitek2s00/skeeper/?locale=en#skeeper'. The user is logged in as 'Laboratory Administrator'. The interface has a top navigation bar with 'Home', 'Settings', 'Help', and 'Quit' buttons. Below this, there are tabs for 'Users' and 'Groups'. The 'Users' tab is active, showing a 'User List' on the left and 'User Details' on the right. The 'User List' contains a scrollable list of users: 'Fse bioMérieux', 'Laboratory Administrator', 'Laboratory Supervisor', and 'Laboratory Technician'. The 'User Details' section has a 'Save' button circled in green. Other buttons include 'Edit', 'New', 'Cancel', 'Delete', and a key icon. The 'User Details' form includes fields for 'First name' (John), 'Last name' (Doe), 'Login name' (DoeJ1), 'E-mail', 'Telephone', 'Preferred language' (English), 'Password' (masked with dots), and 'Confirm password' (masked with dots). There is also a checkbox for 'User account is disabled' and a 'Member of' section with a trash icon and a table with columns 'Group name' and 'Application name'.

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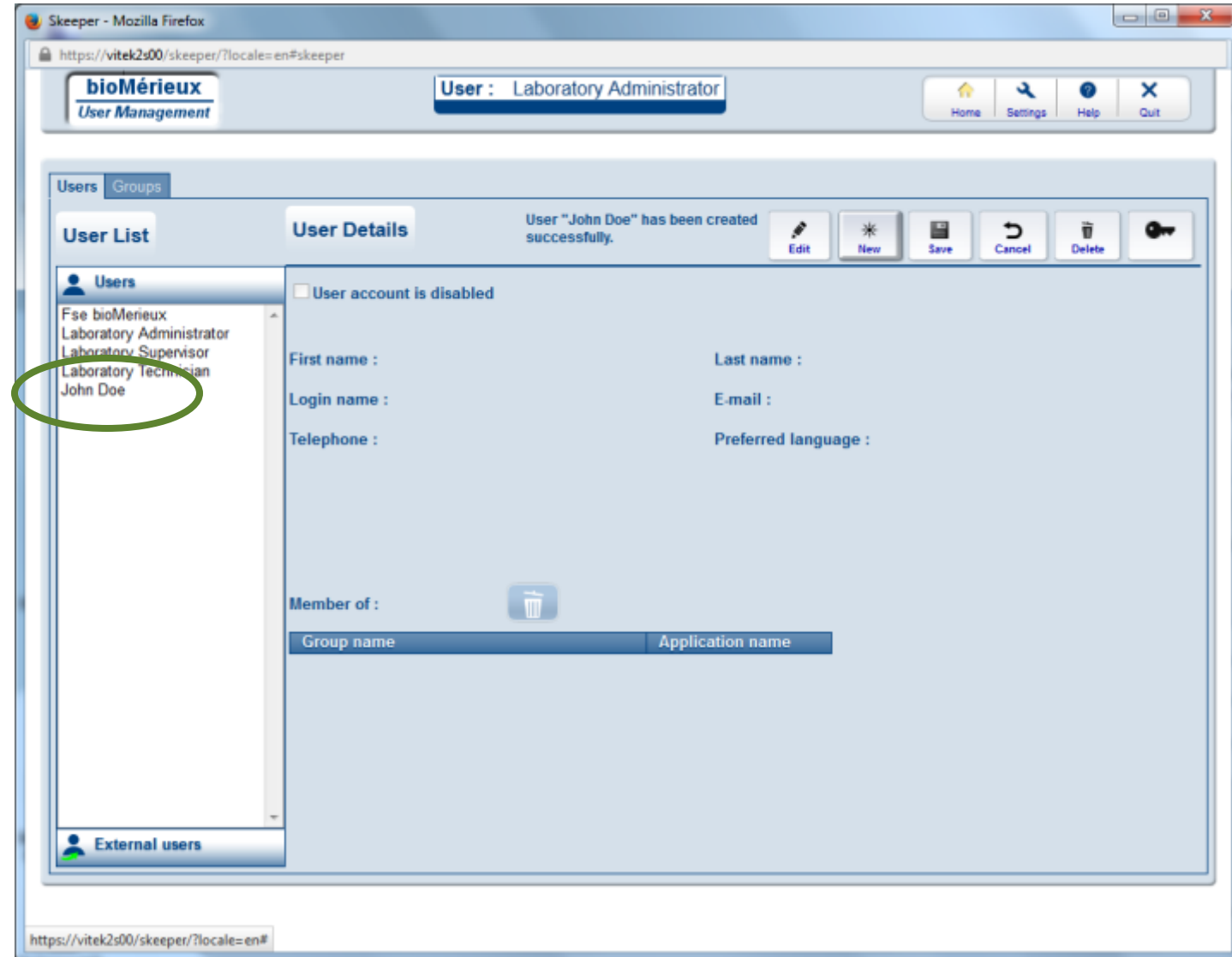The screenshot shows the bioMérieux User Management interface. The browser window is titled 'Skeeper - Mozilla Firefox' and the URL is 'https://vitek2s00/skeeper/?locale=en#skeeper'. The 'User Management' section is active, and the 'User' is set to 'Laboratory Administrator'. The 'Users' tab is selected, showing a 'User List' on the left and 'User Details' on the right. The 'User Details' section includes fields for 'First name', 'Last name', 'Login name', 'E-mail', 'Telephone', and 'Preferred language'. The 'Member of' field is expanded, showing a list of groups. The group 'VITEK2 - Supervisors' is highlighted with a green oval. The 'Add' button is visible at the bottom right.

Group name	Application name
USER MANAGEMENT - UMAAdministrator	
BCICONNECT - Admin	
BCICONNECT - User	
VITEK2 - Supervisors	
VITEK2 - Laboratory Technologists	

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Group name Application name

Select group(s) to add :

USER_MANAGEMENT - UMAAdministrator
BCICONNECT - Admin
BCICONNECT - User
VITEK2 - Supervisors
VITEK2 - Laboratory Technologists

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User List User Details

Users

- Fse bioMérieux
- Laboratory Administrator
- Laboratory Supervisor
- Laboratory Technician
- John Doe

External users

☐ User account is disabled

First name : John Last name : Doe

Login name : DoeJ1 E-mail :

Telephone : Preferred language : English

Member of :

Group name	Application name
Supervisors	VITEK2

Select group(s) to add :

- USER_MANAGEMENT - UMAAdministrator
- BCICONNECT - Admin
- BCICONNECT - User
- VITEK2 - Laboratory Technologists

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